BOOST

THE LIGHT SALON®

ADVANCED LED MASK FACE
MASQUE LED AVANCÉ VISAGE
MASCHERA ADVANCED LED VISO

WELCOME TO THE LIGHT SALON

GETTING STARTED

Please carefully follow the instructions in the user manual for first time use, ensuring you do a patch test (see p.8).

NEED HELP?

For any general queries about your device, controller or warranty.
Email customerservice@thelight-salon.com

THE VIRTUAL LIGHT SALON

Let us treat you and your skin to a one-on-one virtual consultation with our Light Salon Skin Therapists.

Let them help you navigate our signature Boost LED device, recommend additional skincare modalities, deliver advice, and educate on all our offers in-salon.
Email customerservice@thelight-salon.com

www.thelight-salon.com

USING YOUR BOOST LED MASK

**Before Boosting, charge the controller,
using only the supplied cable and plug.**

Do not attach the controller to your LED Mask whilst charging.

Once fully charged (steady green light),
detach the controller from the charging cable.

Your LED Mask will not work if it is attached to the charging cable.

Take your LED Mask and place the longer strap across the top of the head,
shorter strap across the back of the head.

It should feel secure, but comfortable, and not tight.


The LED Mask will naturally sit closer on the top half of the face* and may
lift slightly off the lower half. Shortening the longer strap on the top of your
head will allow the LED Mask to sit higher up.

Attach the controller to your LED Mask.

To turn on your LED Mask firmly press the button on the controller.
Your LED Mask will automatically turn off after 10 minutes.

Store your LED Mask flat, look after it and it will look after you.

*It is normal for a temporary pressure mark to appear on the nose following use.
This will fade up to 30 minutes after treatment.

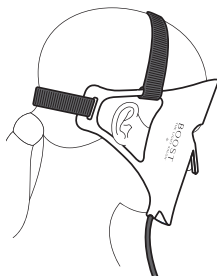BOOST LED MASK REJUVENATING AT-HOME FACIAL

We recommend using your Boost LED Mask 3-5x per week for optimum results, however
it can be used daily according to your skin's needs. Use on its own, on thoroughly
cleansed skin, or follow our guide for a salon experience at-home.

1. CLEANSE	2. SPRITZ	3. REPLENISH
Start with a double cleanse, using our Enzymatic Cleanser & Mask. For a weekly exfoliating mask: Apply a thin layer on dry skin. Leave for 5 to 10 minutes before rinsing thoroughly with water. Pat dry.	Your (optional) second cleanse is a generous spritz of our Cleanse & Recovery Spray.	Apply our Hydrating Peptide Serum*, gently pat into the skin and wait 30 seconds for the serum to absorb. No essential oils, exfoliants or Retinol at this stage.
4. BOOST LED MASK	5. POST CARE	WEEKLY TREAT
Take your LED Mask, secure the back and top straps and deliver a 10 minute dose of light. 	Daytime: *If you prefer, you can apply our Hydrating Peptide Serum post-LED, followed by your SPF. Night-time: Apply your resurfacing or Retinol products, followed by our Hydrating Peptide Serum.	Specifically designed to work with LED; the Hydrogel targets the outer cells whilst the light penetrates the transparent mask to treat the skin from within.
6. REPEAT		
Easy Maintenance, it's that simple. Repeat daily and thank us later.		

CONTENTS

	Welcome	1
	Using your Advanced LED Mask	2
	Rejuvenating At-Home Facial	3
1.0	The Light Salon Boost® Advanced LED Mask	5
	1.1 Intended use	5
2.0	Device Description	5
3.0	Contraindications, Warnings and Precautions	6
4.0	Risks and Benefits	8
	4.1 Potential risks	8
	4.2 Potential benefits	8
5.0	First Time Setup	8
	5.1 Unpacking and charging your Advanced LED Mask	8
	5.2 Connecting the Advanced LED Mask to the controller	8
	5.3 Sensitisation test	9
6.0	Using your Advanced LED Mask device	9
	6.1 Treatment recommendations	9
	6.2 Wearing your Advanced LED Mask	9
	6.3 Carrying out a treatment	10
7.0	Care and Maintenance	10
	7.1 Cleaning your Advanced LED Mask	10
	7.2 Disposal	10
	7.3 Troubleshooting	11
8.0	Labelling	11
9.0	Warranty	12
	9.1 Standard warranty	12
	9.2 Warranty exceptions	12
10.0	Contact Details	13
11.0	Safety and EMC	14
12.0	Appendix	20

1.0 THE LIGHT SALON BOOST® ADVANCED LED MASK

As we age skin becomes thinner, with a drier texture due to dehydration, developing wrinkles due to the loss of collagen and elastic tissue. Aging skin is also much more susceptible to the development of brown spots and pigmentation. The skins dermis can lose from 20-80% of its thickness during the aging process.

The Advanced LED Mask uses two proven wavelengths of light; 633nm red light and 830nm near-infrared light. These wavelengths have been shown in clinical studies to stimulate collagen and elastin production cells and to improve blood flow and tissue oxygenation through a process called photobiomodulation.

The Advanced LED Mask has been designed to deliver light directly to the skin, ensuring that the maximum amount of light hits its cellular targets and is not wasted.

1.1 INTENDED USE

The Advanced LED Mask emits light in the red and near infrared (NIR) region of the light spectrum and is intended for the treatment of cosmetic indications including full-face wrinkles.

2.0 DEVICE DESCRIPTION

The Advanced LED Mask consists of:

1. Silicone flexible LED Mask x1
2. Control unit x1
3. Power supply x1
4. Country specific adaptors x4
5. Long strap x1
6. Short strap x1
7. Protective bag for storing your LED Mask, cable and adaptors x1
8. Eye Shields* x1 Set

* Eye shields clip easily into the eye holes of the LED Mask. Their use is optional as this mask complies to all relevant eye safety standards; they are supplied for comfort only.

The Advanced LED Mask is a home-use wearable LED phototherapy device.

The device consists of a flexible silicone mask that contains light emitting diodes (LEDs) and a controller. The LEDs generate the light. The LED Mask is worn on the face and is held in place by two adjustable Velcro straps.

The controller switches the LEDs ON/OFF and controls power to the LED Mask. The controller contains a rechargeable lithium-ion battery. The controller features two micro-LED indicators to display the battery status:

1. Green Light (Flashing): Device is charging
2. Green Light (Steady): Fully charged
3. White Light: Low battery

The power supply is used to charge the lithium battery and is connected to a suitable electrical outlet via a 2 or 3 pin input socket and wall plug. The power cable is connected to the controller by a standard micro USB-C connector. The Advanced LED Mask will NOT work whilst charging.

This equipment is not intended for diagnostic use or to determine whether treatment is required. The equipment does not require checks on the light output as the LEDs do not dim with age to any practical extent.

The Advanced LED Mask (LED=Light Emitting Diode) emits light energy in the Red and Near-infrared (NIR) region of the light spectrum and is intended to treat facial wrinkles through a non-thermal mechanism called Photobiomodulation.

3.0 CONTRAINDICATIONS, WARNINGS AND PRECAUTIONS

CONTRAINDICATIONS

⚠ DO NOT share the Advanced LED Mask with any other person. The Advanced LED Mask is for single person use only.

• DO NOT use the Advanced LED Mask to treat any other conditions apart from those listed in the indications for use. Advanced LED Mask has not been tested for any other conditions than those listed, and the risk is unknown.

• DO NOT use the Advanced LED Mask if you are pregnant, planning on becoming pregnant or breastfeeding. The Advanced LED Mask has NOT been tested on pregnant or breastfeeding women and therefore the risk to the foetus, neonate or pregnant women is unknown.

• DO NOT use the Advanced LED Mask if you suffer from light induced headaches.

• DO NOT use the Advanced LED Mask if you suffer from any genetic conditions of the eyes. The Advanced LED Mask has NOT been tested on subjects who suffer with genetic conditions of the eyes and therefore the risks are unknown. If you are unsure about any related eye condition and using the Advanced LED Mask contact your ophthalmologist or eye doctor.

• DO NOT use the Advanced LED Mask if you suffer from Lupus erythematosus, photosensitive eczema or Albinism. If you use the Advanced LED Mask to treat Lupus erythematosus, photosensitive eczema or Albinism you may cause a severe skin reaction

• DO NOT use the Advanced LED Mask if you suffer from any photosensitive disorder (sensitisation to light). If you use the Advanced LED Mask and you suffer from a photosensitive disorder you may cause a severe skin reaction.

• DO NOT use the Advanced LED Mask if you are taking any medication that can cause photosensitivity. If you use the Advanced LED Mask and you are taking any medication that can cause photosensitivity you may cause a severe skin reaction.

• For a full list of drugs that can cause photosensitivity, please visit: www.thelight-salon.com/contraindications

⚠ **Photosensitivity** is a common side effect of various medications. These can include certain antibiotics, chemotherapy drugs, and diuretics. If you are unsure about any medication you may be taking consult your healthcare provider.

⚠ **Other substances not listed above can also cause photosensitivity.**

• Common examples of these substances are: St John's Wort, Coal Tar, Deodorants, antibacterial soaps, artificial sweeteners, naphthalene (mothballs), petroleum products, brightening agents found in laundry detergent, and cadmium sulphide (a chemical injected into the skin during tattooing).

WARNINGS

⚠ This device is not intended to treat the eyes or eyelids.

⚠ Before using the Advanced LED Mask you must carry out a sensitisation test (see section 5.3, page 8).

• DO NOT use the Advanced LED Mask when charging.

• DO NOT look directly at the light-emitting diodes (LEDs) when using the Advanced LED Mask.

• DO NOT use if the device or its accessories are damaged, as they may cause injury. Contact customerservice@thelight-salon.com if you suspect your device is damaged.

• DO NOT lie (or place any other heavy weight) on the Advanced LED Mask. The Advanced LED Mask is not designed to withstand such forces. Any damage caused by misuse in this way will invalidate the warranty.

• DO NOT place the device under clothing, blankets, sheets or towels.

• Keep the device away from flammable liquids and aerosol products.

• DO NOT use this product when stacked with other equipment or near interfering equipment, as this may damage the product or cause unexpected treatment results.

• Use only accessories specified or supplied by the manufacturer. Using other accessories may increase electromagnetic emissions or decrease electromagnetic immunity. The device must be used only with the charger and USB cable provided by the manufacturer.

• This equipment has been tested for radiated RF immunity only at specified frequencies. Use near emitters operating at other frequencies may cause unexpected treatment results. (See "Test specifications for ENLOSURE PORT IMMUNITY to RF wireless communications equipment" and "Test specification for enclosure port immunity to proximity magnetic fields.")

• This product is not recommended for use within one month of having facial fillers or similar treatments.

• DO NOT use the Advanced LED Mask if you have experienced sunburn in the past two weeks. Allow your skin to fully recover.

• Keep out of reach of children. This device is not a toy and should not be used or played with by children.

• DO NOT use the Advanced LED Mask if you are sensitive to the Advanced LED Mask.

• DO NOT use the Advanced LED Mask in any way that is not exactly as specified in these instructions for use. Use of the Advanced LED Mask other than exactly as specified in these instructions for use has not been tested and the risk is unknown.

• TURN OFF and REMOVE the Advanced LED Mask immediately if it feels hot or your skin feels uncomfortably warm.

• STOP using the Advanced LED Mask during a treatment if you experience headaches or changes in vision.

• Only use the proprietary hydrogel mask recommended. Other masks are available but have not been tested with the device, and may block light penetration and reduce the effects of the treatment.

• DO NOT use the Advanced LED Mask without first removing creams, lotions or make-up BEFORE treatment. Some ingredients in these products may cause light sensitisation or may reduce the effect of the treatment. If you use Advanced LED Mask without removing these products, you may cause a skin reaction, or you may reduce the effect of the treatment.

• DO NOT use if you have, or suspect you have, cancer. Please consult your doctor.

• NEVER use the Advanced LED Mask near water or where water can be spilt onto the equipment. DO NOT use Advanced LED Mask if you think it has come into contact with water.

• Keep the Advanced LED Mask away from HOT surfaces.

• DO NOT attempt to modify the Advanced LED Mask in any way. Modifying the Advanced LED Mask will void your warranty. Modifying the Advanced LED Mask could cause an electric shock when the system is used again. The system contains no serviceable parts.

GENERAL SAFETY

⚠ The Advanced LED Mask is an electrical device, which under certain circumstances could present an electrical shock hazard to the user. Please read the accompanying user manual thoroughly and follow directions stated in the manual to assure maximum safety during operation.

PRECAUTIONS

⚠ The Advanced LED Mask strap contains a mixture of polyester, nylon and lycra. If you are allergic to any of these materials, you may not be able to use the head strap provided. Speak with your supplier about returning the device.

• Before each use, with the device switched off, all surfaces of the Advanced LED Mask may be cleaned by following the cleaning instructions in Section 7.1, page 9.

• Do NOT use any other cleansing products other than those stated in the cleaning instructions. Use of other cleansing products may damage the Advanced LED Mask.

• Store the Advanced LED Mask at room temperature in dry conditions, away from direct sunlight, ideally in the protective bag.

⚠ DO NOT operate the device near equipment marked with the electromagnetic-frequency warning symbol (e.g. cellular towers).

⚠ DO NOT use the device in environments where other devices emit unshielded electromagnetic energy. Medical electrical equipment can be affected by portable or mobile RF communication equipment.

⚠ DO NOT expose the device to MRI, diathermy, or electrocautery. The device is MR-unsafe. Do not bring it into MR (magnetic resonance) environments.

4.0 RISKS AND BENEFITS

4.1 POTENTIAL RISKS

Immediately after treatment skin may appear red. Skin should return to normal within 1-2 hours.

Temporary headaches immediately following the light treatment. Headaches can be treated with an appropriate over the counter analgesic and rehydration.

4.2 POTENTIAL BENEFITS

• People who have undergone a course of 12 Advanced LED Mask treatments report a reduction in fine lines and periorbital wrinkles, and improvements in texture.

5.0 FIRST TIME SETUP

5.1 UNPACKING AND CHARGING YOUR ADVANCED LED MASK

The Advanced LED Mask uses a controller which contains a lithium-ion battery. When the controller is plugged in and charging; do not attempt to use the Advanced LED Mask as it will not illuminate.

The lifetime of the lithium battery can be maximised by allowing it to fully charge and then fully discharge during regular use.

It will take approximately 4 hours to fully charge the controller.

• Unpack the Advanced LED Mask and check off the parts described in Section 2, page 5.

• Before use, fully charge the battery.

• To charge the Advanced LED Mask battery insert the USB-C cable into the controller.

• Insert the power supply into a suitable electrical outlet and switch on.

• The micro-LEDs in the controller will start to flash to indicate charging. The lithium-ion battery is fully charged when the micro-LED indicator displays a steady green light.

5.2 CONNECTING THE ADVANCED LED MASK TO THE CONTROLLER

Connect the Advanced LED Mask by inserting the USB connector into the USB port of the controller.

5.3 SENSITISATION TEST

Before carrying out your first treatment you must carry out a sensitisation test. Follow the step by step guide below:

1. Connect the Advanced LED Mask to the controller.
2. Place the Advanced LED Mask over your inner arm.
3. Switch ON the controller and allow the Advanced LED Mask to run for its full 10-minute treatment program.
4. At any time if your skin feels HOT or uncomfortable STOP using the device immediately.
5. After the treatment wait at least 6 hours. If your skin on your inner arm is red or itchy, the Advanced LED Mask is unsuitable for you. Return the device to the retailer within (14) days of purchase.

If after 6 hours there is no reaction please repeat the test on your other forearm, however this time we want to ensure that your routine skin cleanser is compatible with light.

On the rare occasion a cleanser can cause photosensitivity due to specific combinations of active ingredients so testing it with the light means that you can be assured of compatibility.

6. Cleanse your arm with your cleanser of choice and remove product thoroughly with warm water.
7. Place the Advanced LED Mask to the cleansed area.
8. Switch ON the controller and allow the Advanced LED Mask to run for its full 10-minute treatment program.
9. At any time if your skin feels HOT or uncomfortable STOP using the device immediately.
10. After the treatment wait at least 6 hours. If your skin on your inner arm is red or itchy, email customerservice@thelight-salon.com with your cleanser and ingredient list.

6.0 USING THE ADVANCED LED MASK

6.1 TREATMENT RECOMMENDATIONS

• For best results treat for 3-5 days a week for a minimum of 6 weeks.

• Maintenance treatments should be carried out as required, according to lifestyle and skin condition.

6.2 WEARING THE ADVANCED LED MASK



• The Advanced LED Mask comes with two adjustable Velcro straps that allow you to comfortably wear the LED Mask.

• To fit the LED Mask, thread the shorter strap through the inlet holes on either side of the LED Mask as shown in the images below.

• The longer strap fits through the top of the loops on either side of the LED Mask, as shown in the images, and sits on the top of your head.

• Loop the straps around your head and secure firmly.

• Adjust if necessary.

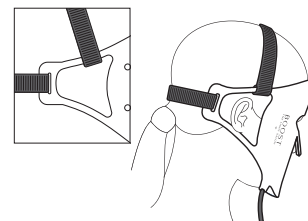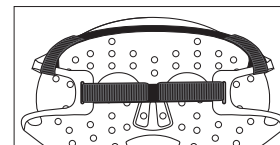6.3 CARRYING OUT A TREATMENT

1. Cleanse and gently exfoliate your skin to remove any traces of make-up or skincare.
2. Place your Advanced LED Mask in position, and secure straps around the head (longer strap sits on top of head, shorter strap secured around the back of head).
3. Connect the Advanced LED Mask to the controller and switch on the controller.

▲ TOP TIP: Tuck the controller into the shoulder strap of your clothing, or place in a pocket if choosing to wear whilst standing in order to go about your day.

4. Wear for 10 minutes (automatic shut-off). You may stop the treatment at any time by pressing the ON/OFF button on the controller.
5. Remove your Advanced LED Mask and finish with skincare products of your choice.
6. Clean the LED Mask using a water-based wipe or a lint-free damp clean cloth. Ensure that the device is NOT being charged when cleaning.
7. Repeat 3-5 times a week for 6 weeks.

For more information on recommended products, go to www.thelight-salon.com

7.0 CARE AND MAINTENANCE

- The lifetime of the lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the Advanced LED Mask at room temperature in dry conditions. Use the protective bag provided.
- Avoid storing the Advanced LED Mask in dusty conditions, in direct sunlight or in damp and humid conditions.
- The Advanced LED Mask contains no serviceable parts.
- The Advanced LED Mask should always be stored FLAT. The flexible technology is designed for flexible comfort on the face, not for storage.
- Before each use, with the device switched off, clean the pads following the cleaning instructions below.

7.1 CLEANING THE ADVANCED LED MASK

DO NOT use any other cleansing products as these may damage the plastics.

- Clean the transparent (face side) of the LED Mask with a water-based wipe or damp clean cloth.
- Clean the rear of the LED Mask, and the entire surface of the controller with a water-based wipe or damp clean cloth.
- Place the LED Mask on a clean dry surface to air dry. Visually inspect the surfaces to ensure they are clean. Stubborn stains can be removed by wiping gently with 70% alcohol.
- If you follow the care and maintenance for your Advanced LED Mask, the Advanced LED Mask has a minimum lifetime of 2 years

7.2 DISPOSAL

- ▲** DO NOT dispose of any part of the Advanced LED Mask at a domestic waste facility.

 The Advanced LED Mask controller contains a lithium-ion battery. For protection of the environment, please package and return this unit to the manufacturer or their representative at the end of its working life or safely dispose of the device following your local guidelines.

7.3 TROUBLESHOOTING

Problem	Possible Cause	Action
No red or infrared light from the Advanced LED Mask	Not charged	Plug in the controller to the charger making sure that the connectors are attached and leave charging until the three indicator LEDs on the controller remain lit and do NOT flash.
No red or infrared light from the Advanced LED Mask	Device connected to the mains	Unplug micro USB connector from controller
LEDs do NOT turn OFF after 10 minutes	Device is faulty	Return to your supplier
LEDs flicker during treatment	Device is faulty	Return to your supplier

▲ If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please return to your provider.

▲ Do NOT attempt to investigate any further.

▲ Do NOT attempt to modify the device in any way. The Advanced LED Mask has NO parts you can repair yourself.

TECHNICAL SPECIFICATION

- Program time: 10 minutes fixed
- Operating Temperature: 0°C to 40°C
- Storage and Transport Temperature: 20°C to 60°C
- Operating Humidity: 30% to 80% without condensation
- Storage and Transport Humidity: 30% to 90% without condensation
- Electrical Input from power adaptor: 100-240V-50/60Hz 0.4A
- Electrical Output: 5V, 1.0A
- 3.6V, 9.18Wh, 2550mAh

8.0 LABELLING

The labelling supplied with the device contains the following information:

- Date of manufacture
- Serial number and batch code
- Operating environment

RECYCLING

 This symbol means that this product shall not be disposed of with normal household waste (2012/19/EU).

We advise you to take your product to an official collection point to have a professional remove the rechargeable battery. Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent damage to the environment and human health.

EXPLANATION OF SYMBOLS

There are several safety warning labels on the Advanced LED Mask and associated packaging and in the user manual. Not following these warnings may result in one or more of the following hazards: electrical shock and/or fire.

-  This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
-  This symbol on the device label means: Read the user manual before you start using the device.
-  This symbol means: that the Advanced LED Mask is specified as applied part type BF, B indicates that it is used on the Body and F indicates that it is not connected to the mains voltage.
-  Indicates the device is CE Marked and authorised to sell in Europe (Low Voltage Directive and Electromagnetic Compatibility Directive).
-  This symbol implies conformity to UK Directives on health, safety, and environmental protection standards. UKCA stands for UK Conformity Assessed.
-  This symbol means: 'Manufactured by' and provides the name and address of the legal manufacturer.
-  This symbol indicates input for the power supply. This symbol means that the device may only be charged with the adapter that is supplied in the package. Usage of other adapters may not be safe and could lead to electrical shock.
-  This symbol means: Keep dry.
-  This symbol indicates the serial number of the device.
-  Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission, and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
-  This symbol refers to the batch code of the Advanced LED Mask.
-  This symbol indicates the USB connector.

-  This symbol means that the device uses a lithium-ion battery.
-  This symbol means: Made in China.
-  This symbol means: Keep the device away from direct sunlight.
-  This symbol means the upper and lower operating temperatures that the Advanced LED Mask can be safely used at.
-  This symbol means the upper and lower operating humidity that the Advanced LED Mask can be safely used at.
-  This symbol means that the cardboard packaging that the Advanced LED Mask came in can be recycled.

IP22 This symbol means that the device is protected against a solid object greater than 12.5mm such as a finger and against vertically falling drops of water at up to 15 degrees from vertical.

9.0 WARRANTY

Please register your warranty at www.thelight-salon.com

9.1 STANDARD WARRANTY

Under normal use this product is guaranteed for 2 years after its original purchase date against manufacturing defects in materials or workmanship.

9.2 WARRANTY EXCEPTIONS

The lithium battery is guaranteed for 1 year after its original purchase date against manufacturing defects in materials or workmanship.

This warranty extends only to the original purchaser.

This warranty does not cover damage caused by the following:

- Accident, misuse, abuse or alteration
- Servicing by unauthorized persons
- Connection to incorrect power sources
- Use with unauthorized accessories

WHAT TO DO IF YOUR PRODUCT IS FAULTY

If the product should become faulty within the warranty period, contact your place of purchase.

The Light Salon Boost Limited reserves the right to replace a defective product with the most compatible product currently available.

The Light Salon Boost Limited shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to one year from date of original purchase.

This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you.

This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed.

-  Users of this device can obtain benefits if all information in the user manual is read and comprehended. All the instructions in the user manual must be followed.

Please contact us on our Customer Services email customerservice@thelight-salon.com if you have any comments about THE LIGHT SALON BOOST® Advanced LED Mask or User Guide.

10.0 CONTACT DETAILS

Registered office
71-75 Shelton Street, Covent Garden,
London, WC2H 9JQ, United Kingdom

www.thelight-salon.com

 The information in this document is proprietary to The Light Salon Limited, and supplied on the understanding that it shall be kept confidential to recipients and that it shall not be reproduced and/or used for commercial purposes without the written consent of THE LIGHT SALON BOOST® LIMITED.

 **Shenzhen Kaiyan Medical Equipment Co. Ltd**
Building*3 and Building*5
40th of Fuxin Street, Huaide Community
Fuyong Town, Baoan District, Shenzhen,
Guangdong 518103, China

EU Authorised Representative: Light Tree Ventures Europe B.V.
Laan van Ypenburg 108, 2497 GC,
The Hague, The Netherlands
regulatory@lighttreeventures.com

UK Authorised Representative: VH Consulting Limited
Office 147, Centurion House, London Road,
Staines-upon-Thames, Surrey, TW18 4AX
United Kingdom
h2yhuk@gmail.com

11.0 SAFETY AND EMC

This device is internally powered equipment with type BF applied part. It complies with Medical Electrical Safety Standards (IEC 60601-1, IEC 60601-1-11, IEC 60601-2-57).

This device also complies with Medical EMC Standard (IEC 60601-1-2).

The device's essential performance is: The main functions of the subject's device, such as button function, indicator light display, LED light at the application site, charging and overvoltage charging protection, can be maintained in normal operation.

If the functions/performance are lost or degraded due to electromagnetic (EM) disturbances, the user may notice issues such as the controller button becoming unresponsive, indicator lights and LED lights flickering, becoming non-functional, or showing a sudden increase in power.*

*Please refer to troubleshooting (7.3) to resolve any such issue. If the issue cannot be resolved, it may be caused by electromagnetic (EM) disturbances. Keep the device away from the source of EM emitters and ensure it operates in a recommended electromagnetic environment.

GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC EMISSIONS		
The ABC is intended for use in the electromagnetic environment specified below. The customer or the user of the ABC should ensure that it is used in such an environment.		
Emissions	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The ABC is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	$P_d \leq 1.0$ $P_n \leq 0.65$ $d_{(0.1\text{ Hz})} \leq 0.5$ $d_{\text{max}} \leq 4\%$ $d_r \leq 3.3\%$	

GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY			
The ABC is intended for use in the electromagnetic environment specified below. The customer or the user of the ABC should ensure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV 100 kHz repetition frequency	± 2 kV 100 kHz repetition frequency	Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT.
Surge IEC 61000-4-5	± 1 kV Line-to-line	± 1 kV Line-to-line	Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°) 0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°) 0 % UT; 250/300 cycles	0 % UT; 0.5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°) 0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°) 0 % UT; 250/300 cycles	Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT. If the user of the ABC requires continued operation during power mains interruptions, it is recommended that the ABC be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a HOME HEALTHCARE ENVIRONMENT.
NOTE: UT is the a.c. mains voltage prior to the application of the test level			

GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY

The ABC is intended for use in the electromagnetic environment specified below.
The customer or the user of the ABC should ensure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Conducted RF IEC 61000-4-6	3 Vrms 0.15 MHz to 80 MHz 6 Vrms in ISM and amateur radio bands between 0.15 MHz to 80 MHz	3 Vrms 0.15 MHz to 80 MHz 6 Vrms in ISM and amateur radio bands between 0.15 MHz to 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of the ABC, including cables, than the recommended separation distance calculated from the equation applicable to the transmitter's frequency. $d = 1,2\sqrt{P}$
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz to 2.7 GHz 80 % AM at 1 kHz	$d = 1,2\sqrt{P}$ (80 MHz to 800 MHz) $d = 2,3\sqrt{P}$ (800 MHz to 2.7 GHz) Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer, and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b Interference may occur in the vicinity of equipment marked with the following symbol: 

Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a) Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast, and TV broadcast, cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength at the location where the ABC is used exceeds the applicable RF compliance level above, the ABC should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the ABC.

b) Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

RECOMMENDED SEPARATION DISTANCES BETWEEN PORTABLE AND MOBILE RF COMMUNICATIONS EQUIPMENT AND THE ABC

The ABC is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the ABC can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the ABC as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1,2\sqrt{P}$	80 MHz to 800 MHz $d = 1,2\sqrt{P}$	800 MHz to 2.7 GHz $d = 2,3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

TEST SPECIFICATIONS FOR ENCLOSURE PORT IMMUNITY TO RF WIRELESS COMMUNICATIONS EQUIPMENT						
Test Freq. (MHz)	Band (MHz)	Service	Modulation	Max. Power (W)	Distance (m)	Immunity Test Level (V/m)
385	380-390	TETRA 400	Pulse modulation 18 Hz	1.8	0.3	27
450	430-470	GMRS 460, FRS 460	FM ± 5 kHz deviation 1 kHz sine	2	0.3	28
710	704-787	LTE 13, 17 Band	Pulse modulation 217 Hz	0.2	0.3	9
745						
780						
810	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE 5 Band	Pulse modulation 18 Hz	2	0.3	28
870						
930						
1720	1700-1990	GSM 1800, CDMA 1900, GSM 1900, DECT, LTE 1, 3, 4, 25 Band, UMTS	Pulse modulation 217 Hz	2	0.3	28
1845						
1970						
2450	2400-2570	Bluetooth, WLAN 802.11 b/g/n, RFID 2450, LTE 7 Band	Pulse modulation 217 Hz	2	0.3	28
5240	5100-5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9
5500						
5785						

TEST SPECIFICATION FOR ENCLOSURE PORT IMMUNITY TO PROXIMITY MAGNETIC FIELDS

Test Frequency	Modulation	Immunity Test Level (A/m)
30 kHz	CW	8
134.2 kHz	Pulse modulation 20.1 kHz	65
13.56 MHz	Pulse modulation 50 kHz	7.5
a) This test is applicable only to ME EQUIPMENT and ME SYSTEMS intended for use in the HOME HEALTHCARE ENVIRONMENT. b) The carrier shall be modulated using a 50 % duty cycle square-wave signal. c) r.m.s., before modulation is applied.		

12.0 APPENDIX

According to MDD 93/42/EEC, the Advanced LED Mask is classified as a cosmetic device within the European union. However, for purposes of safety we have conducted additional testing so that you can be ensured of the safety and efficacy of the Advanced LED Mask.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The Advanced LED Mask has been tested to or complies with the following standards:

IEC 60601-1:2005+AMD1:2012+AMD2:2020 Medical electrical equipment. General requirements for basic safety and essential performance.

IEC 60601-1-2:2014+A1:2020 Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Collateral Standard: Electromagnetic disturbances - Requirements and tests.

IEC 60601-1-6:2010+AMD2:2020 Medical electrical equipment - Part 1-6: General requirements for basic safety and essential performance - Collateral standard: Usability.

IEC 60601-1-11:2015+AMD1:2020 Medical electrical equipment - Part 1-11: General requirements for basic safety and essential performance - Collateral Standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment.

IEC 60601-2-57:2023 Medical Electrical Equipment - Part 2-57: Particular requirements for the basic safety and essential performance of non-laser light source equipment intended for therapeutic, diagnostic, monitoring, cosmetic and aesthetic use.

IEC 62471:2016 Photobiological Safety of Lamps and Lamp Systems.

IEC 62366-1:2015 Medical devices - Part 1: Application of usability engineering to medical devices

ISO 10993-1:2018 Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process.
ISO 15223-1:2021 Medical Devices - Symbols to be used with information to be supplied by the manufacturer - General requirements.

EN ISO 14971:2019 Medical devices - Application of risk management to medical devices.

IEC 62133-2:2017+AMD1:2021 Secondary cells and batteries containing alkaline or other non-acid electrolytes - Safety requirements for portable sealed secondary lithium cells, and for batteries made from them, for use in portable applications - Part 2: Lithium systems.

BS EN 55014-1:2006+A2:2011 Electromagnetic compatibility. Requirements for household appliances, electric tools and similar apparatus - Part 1: Emission.

BS EN IEC 55014-2:2021 or EN IEC 55014-2:2021 Electromagnetic compatibility - Requirements for household appliances, electric tools and similar apparatus - Part 2: Immunity - Product family standard.

RoHS Directive 2011/65/EU Annex II amending Directive (EU)2015/863 The Restriction of the use of certain hazardous substances in electrical and electronic equipment.

The Advanced LED Mask is cleared for marketing in the United States of America K221775.

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